



**IFN-gamma (pig)**

A brand name of



**For laboratory research use only. Not for human or veterinary diagnostic use.**

**Bertin Pharma also markets pre-analytical products, EIA kits, antibodies & biochemicals for:**

- ▶ **Cardiology / Hypertension**
- ▶ **Diabetes / Obesity**
- ▶ **Endocrinology / Metabolism**
- ▶ **Inflammation**
- ▶ **Pharmacology**
- ▶ **Psychopharmacology**
- ▶ **Nitric Oxide**
- ▶ **Oncology / Apoptosis**
- ▶ **Oxidative injury**
- ▶ **Cell signaling**
- ▶ **Drug metabolism**

**Do not hesitate to contact our after-sales services for further information at [bioreagent@bertinpharma.com](mailto:bioreagent@bertinpharma.com)**



**IFN-gamma (pig)  
Enzyme Immunoassay kit  
#A05413.96 wells**

For research laboratory use only  
Not for human diagnostic use

This assay has been developed & validated  
by Bertin Pharma



Fabriqué en France  
Made in France

#A11413  
Version: 0117

## Table of contents

▶ Precaution for use	6
▶ Background	7
▶ Principle of the assay	8
▶ Materials and equipment required	10
▶ Sample collection and preparation	11
▶ Reagent preparation	12
▶ Assay procedure	15
▶ Data analysis	20
▶ Acceptable range	21
▶ Typical results	22
▶ Characteristics	24
▶ Troubleshooting	25
▶ Bibliography	26

**96 wells**  
**Storage: +4°C**  
**Expiry date: stated on the package**

This kit contains:

Designation	Colour of cap	Item #	Quantity per kit	Form
IFN-gamma precoated 96-well Strip Plate	Blister with zip	A08413.1 ea	1	-
Streptavidin Poly_HRP Tracer	Green	A04410.100 dtn	1	Liquid
IFN-gamma (pig) Biotin-labelled Antibody	Red	A03413.100 dtn	1	Liquid
IFN-gamma (pig) Standard	Blue with red septum	A06413.1 ea	2	Lyophilised
Poly_HRP EIA Buffer	Grey / Blue	A07410.1 ea	1	Lyophilised
Wash Buffer	Silver	A17000.1 ea	1	Liquid
Tween 20	Transparent	A12000.1 ea	1	Liquid
HRP Substrate Solution	Black	A09034.100 dtn	1	Liquid
HRP Stop Solution	Yellow	A22410.100 dtn	1	Liquid
Instruction Booklet	-	A11413.1 ea	1	-
Well cover Sheet	-	-	1	-

Each kit contains sufficient reagents for 96 wells. This allows for the construction of one standard curve in duplicate and the assay of 37 samples in duplicate.

## ► **Precaution for use**

**Users are recommended to carefully read all instructions for use before starting work.**

Each time a new pipette tip is used, aspirate a sample or reagent and expel it back into the same vessel. Repeat this operation two or three times before distribution in order to equilibrate the pipette tip.

- For research laboratory use only
- Not for human diagnostic use
- Do not pipet liquids by mouth
- Do not use kit components beyond the expiration date
- Do not eat, drink or smoke in area in which kit reagents are handled
- Avoid splashing

HRP Stop Solution and HRP Substrate Solution are harmful solutions. To avoid any contact, wear eye, hand, face and clothing protection when handling these.

Wearing gloves, laboratory coat and glasses is recommended when assaying kit materials and samples.

## ► **Temperature**

Unless otherwise specified, all the experiments are done at room temperature (RT), that is around +20°C. Working at +25°C or more affects the assay.

## ► Background

### ▷ IFN-gamma (pig)

Also known as type II interferons, IFN-gamma is a glycoprotein of 146 amino acids [1].

IFN-gamma is a cytokine critical for the innate and adaptive immunity. It is produced predominantly by natural killer (NK) and natural killer T (NKT) cells as part of the innate immune response and by T helper (Th) cells and cytotoxic T lymphocytes (CTLs) as part of the adaptive immunity [2, 3, 4].

IFN-gamma has antiviral, immunoregulatory and anti-tumor properties [5].

IFN-gamma:

- > promotes the NK cell activity,
- > increases the antigen presentation by action on the lysosome activity of macrophages,
- > promotes Th1 differentiation, and suppresses Th2 differentiation which would cause a humoral (antibody) response,
- > modulates the production of IgG2a and IgG3 from activated plasma B cells,
- > increases expression of class I MHC molecules as well as class II MHC on antigen-presenting cells,
- > promotes adhesion and binding required for leukocyte migration.

IFN-gamma interacts with other cytokines, either in a synergistic (e.g., TNF) or antagonistic (e.g., IL-4) fashion.

## ▶ Principle of the assay

This Enzyme Immunometric Assay (EIA/ELISA) is based on a sandwich technique. The wells of the plate supplied are coated with a polyclonal antibody specific to IFN-gamma (pig).

IFN-gamma (pig) introduced into the wells (standard or sample) will be bound by the polyclonal antibody coated on the plate and is then detected by a second polyclonal antibody tagged with biotin also specific for IFN-gamma (pig).

The two antibodies then form a sandwich by binding on different parts of IFN-gamma (pig).

The sandwich is immobilised on the plate so reagents in excess may be washed away.


The immunological complex is revealed by the interaction between biotin and streptavidin labelled with HRP (Tracer).

The concentration of IFN-gamma (pig) is determined by measuring the enzymatic activity of immobilized Tracer using TMB. The Tracer acts on TMB to form a yellow compound after the reaction has been stopped.

The intensity of the colour, which is determined by spectrophotometry at 450 nm, is proportional to the amount of IFN-gamma (pig) present in the well during the immunological incubation.



The principle of the assay is summarised below:

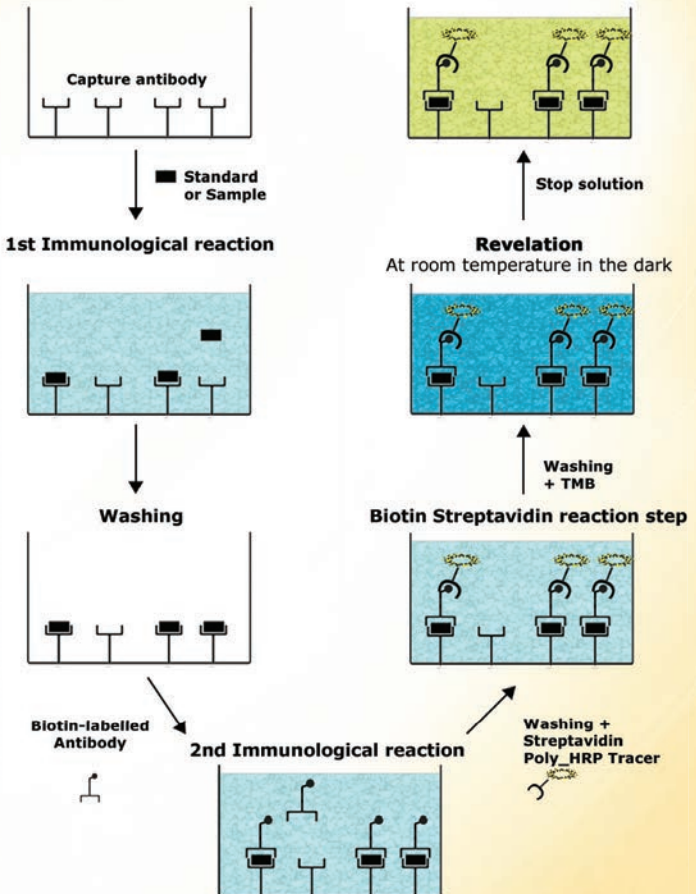

## ► **Materials and equipment required**

In addition to standard laboratory equipment, the following material is required:

For the assay:

- > Precision micropipettes (20 to 1000  $\mu$ L)
- > Spectrophotometer plate reader (450 nm filter)
- > Microplate washer (or washbottles)
- > Orbital microplate shaker
- > Multichannel pipette and disposable tips 30-300 $\mu$ L
- > UltraPure water #A07001.1L
- > Polypropylene tubes



Water used to prepare all ELISA reagents and buffers must be UltraPure, deionized & free from organic contaminant traces.

Do not use distilled water, HPLC-grade water or sterile water.

- > UltraPure water may be purchased from Bertin Pharma (item #A07001.1L).

## ▶ **Sample collection and preparation**

This assay may be used to measure IFN-gamma (pig).



It is the responsibility of the user to check the compatibility of the assay with the study matrix.

## ▷ **General precautions**

- All samples must be free from organic solvents prior to assay.
- Samples should be assayed immediately after collection or should be stored at -20°C or at -80°C.

## ▶ **Reagent preparation**

All reagents need to be brought to room temperature (around +20°C) prior to the assay.

### ▶ **Poly\_HRP EIA Buffer**

Reconstitute Poly\_HRP EIA Buffer #A07410 with 25 mL of UltraPure water. Allow it to stand 5 minutes until completely dissolved and then mix thoroughly by gentle inversion.

*Stability at 4°C: 1 month.*



Before use, filter the Buffer on 0,22 µm filter.

### ▶ **IFN-gamma (pig) Standard**

Reconstitute the IFN-gamma (pig) Standard vial #A06413 with 1 mL of UltraPure water. Allow it to stand 5 minutes until completely dissolved and then mix thoroughly by gentle inversion.

The concentration of the first standard (S1) is 10 ng/mL. Prepare seven polypropylene tubes (for the seven other standards) and add 500 µL of Poly\_HRP EIA Buffer into each tube. Then prepare the standard concentrations by serial dilutions as follow:

Standard	Volume of Standard	Volume of Poly_HRP EIA Buffer	Standard concentration
S1	-	-	10.00 ng/mL
S2	500 µL of S1	500 µL	5.00 ng/mL
S3	500 µL of S2	500 µL	2.50 ng/mL
S4	500 µL of S3	500 µL	1.25 ng/mL
S5	500 µL of S4	500 µL	0.62 ng/mL
S6	500 µL of S5	500 µL	0.31 ng/mL
S7	500 µL of S6	500 µL	0.16 ng/mL
S8	500 µL of S7	500 µL	0.08 ng/mL

*Stability at 4°C: within the day.*

### ▷ IFN-gamma (pig) Biotin-labelled Antibody

The IFN-gamma Biotin-labelled Antibody is supplied concentrated 10 times. Calculate the volume needed (number of wells multiplied by 0.1 mL) and then dilute the IFN-gamma Biotin-labelled Antibody solution #A03413 with the appropriate volume of Poly\_HRP EIA Buffer.

Example: for 40 wells you need 4 mL of IFN-gamma biotin-labelled Antibody (40 x 0.1 mL), add 0.4 mL of IFN-gamma Biotin-labelled Antibody in 3.6 mL of Poly\_HRP EIA Buffer.

*Stability of diluted antibody at +4°C: within the day.*

## **Wash Buffer**

Dilute 2 mL of concentrated Wash Buffer #A17000 with 800 mL of UltraPure water. Add 400  $\mu$ L of Tween 20 #A12000.

Use a magnetic stirring bar to mix the content.

*Stability at +4°C: 1 month.*

## ▶ **Assay procedure**

It is recommended to perform the assays in duplicate following the instructions hereafter.

### ▶ **Plate preparation**

Prepare the Wash Buffer as indicated in the reagent preparation section.

Open the plate packet and select the sufficient strips for your assay. Put unused strips back in the zip lock bag with the dessicant pocket and properly close it.

*Stability at +4°C: 1 month.*

Rinse each well 5 times with Wash Buffer (300 µL/well).

Just before distributing the reagents and samples, remove the buffer from the wells by inverting the plate and shaking out the last drops on a paper towel.

### ▶ **Plate set-up**

A plate set-up is suggested on the following page.

The contents of each well may be recorded on the template sheet provided at the end of this technical booklet.

### ▶ **Pipetting the reagents**

All samples and reagents must reach room temperature prior to performing the assay.

Use different tips to pipet all the reagents.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Bk	S7	S3	*	*	*	*	*	*	*	*	*
B	Bk	S7	S3	*	*	*	*	*	*	*	*	*
C	Bk	S6	S2	*	*	*	*	*	*	*	*	*
D	NSB	S6	S2	*	*	*	*	*	*	*	*	*
E	NSB	S5	S1	*	*	*	*	*	*	*	*	*
F	NSB	S5	S1	*	*	*	*	*	*	*	*	*
G	S8	S4	*	*	*	*	*	*	*	*	*	*
H	S8	S4	*	*	*	*	*	*	*	*	*	*

Bk : Blank

NSB : Non Specific Binding

S1-S8 : Standards 1-8

\* : Samples

Before pipetting, equilibrate the pipette tips in each reagent.  
Do not touch the liquid already in the well when expelling with the pipette tip.

> **Poly\_HRP EIA Buffer**

Dispense 100  $\mu$ L to NSB wells.

> **IFN-gamma (pig) Standard/Sample**

Dispense 100  $\mu$ L of each of the eight standards (S8 to S1) in duplicate to appropriate wells.



Start with the lowest concentration standard (S8) and equilibrate the tip in the next higher standard before pipetting.

### ▷ **Incubating the plate**

Cover the plate with the cover sheet and incubate 120 minutes at room temperature, shaking at 300 rpm.

### ▷ **Washing the plate**

Rinse each well 5 times with the Wash Buffer (300  $\mu$ L/well). Just before distributing reagents, remove the buffer from the wells by inverting the plate and shaking out the last drops on a paper towel.

### ▷ **Pipetting the reagents**

#### ➤ **IFN-gamma (pig) Biotin-labelled antibody**

Dispense 100  $\mu$ L to each well, except Blank (Bk) wells.

### ▷ **Incubating the plate**

Cover the plate with the cover sheet and incubate 60 minutes at room temperature, shaking at 300 rpm.

### ▷ **Washing the plate**

Rinse each well 5 times with the Wash Buffer (300  $\mu$ L/well). Just before distributing reagents, remove the buffer from the wells by inverting the plate and shaking out the last drops on a paper towel.

## ▷ **Pipetting the reagents**

### > **Streptavidin Poly\_HRP Tracer**

Dispense 100  $\mu$ L to each well, except Blank (Bk) wells.

## ▷ **Incubating the plate**

Cover the plate with the cover sheet and incubate 30 minutes at room temperature, shaking at 300 rpm.

## ▷ **Developing and reading the plate**

- > Rinse each well 5 times with 300  $\mu$ L of Wash Buffer. At the end of the last washing step, empty the plate and blot the plate on a paper towel to discard any trace of liquid.
- > Add 100 $\mu$ L of HRP Substrate Solution to each well. Incubate the plate in the dark at room temperature without shaking. For the time, look at the lot specific Quality Control Sheet (QCS). In general, revelation time is 10 min.
- > Add 100 $\mu$ L of HRP Stop Solution to each well.
- > Wipe the bottom of the plate with a paper towel, and make sure that no liquid has splashed outside the wells.
- > Read the plate at 450 nm (yellow color).

Enzyme Immunoassay Protocol (volumes are in $\mu\text{L}$ )				
	Blank	NSB	Standard	Sample
Poly_HRP EIA Buffer	-	100	-	-
Standard	-	-	100	-
Sample	-	-	-	100
Cover plate, incubate 120 minutes at room temperature under orbital shaking at 300 rpm				
Wash strips 5 times with 300 $\mu\text{L}$ /well Discard liquid from the wells & dry on absorbent paper				
Biotin-labelled Antibody	-	100		
Cover plate, incubate 60 minutes at room temperature under orbital shaking at 300 rpm				
Wash strips 5 times with 300 $\mu\text{L}$ /well Discard liquid from the wells & dry on absorbent paper				
Streptavidin Poly_HRP Tracer	-	100		
Cover plate, incubate 30 minutes at room temperature under orbital shaking at 300 rpm				
Wash strips 5 times with 300 $\mu\text{L}$ /well Discard liquid from the wells & dry on absorbent paper				
HRP Substrate Solution	100			
Incubate the plate in the dark without agitation				
HRP Stop Solution	100			
Read the plate at 450 nm				

## ► Data analysis

Make sure that your plate reader has subtracted the absorbance readings of the blank wells (absorbance of HRP Substrate Solution alone) from the absorbance readings of the rest of the plate. If it is not the case, please do it.

- Calculate the average absorbance for each NSB, standard and sample.
- For each standard, plot the absorbance on y axis versus the concentration on x axis. Draw a best-fit line through the points.
- To determine the concentration of your samples, find the absorbance value of each sample on the y axis.
- Read the corresponding value on the x axis which is the concentration of your unknown sample.
- Samples with a concentration greater than 10 ng/mL should be re-assayed after dilution in Poly\_HRP EIA Buffer.
- Most plate readers are supplied with a curve-fitting software capable of graphing these data (4-parameter logistic fit 4PL). If you have this type of software, we recommend using it. Refer to it for further information.

## ► **Acceptable range**

- > NSB absorbance  $\leq 0.060$  A.U.
- > Limit of detection  $\leq 0.08$  ng/mL


## ► Typical results

The following data are for demonstration purpose only. Your data may be different and still correct.

These data were obtained using all reagents as supplied in this kit under the following conditions: 10 minutes developing at room temperature, reading at 450 nm. A 4 parameter logistic fitting was used to determine the concentrations.

Standard	IFN-gamma (pig) ng/mL	Absorbance A.U.
S1	10.00	1.513
S2	5.00	0.911
S3	2.50	0.516
S4	1.25	0.303
S5	0.62	0.172
S6	0.31	0.108
S7	0.16	0.070
S8	0.08	0.060
NSB	0.00	0.041

Typical IFN-gamma (pig) standard curve

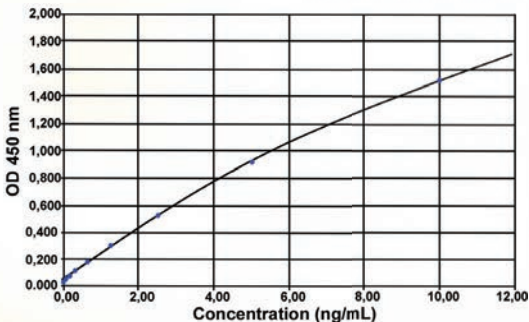

## ► Characteristics

### > Cross-reactivity

Species	Cross-reactivity
Recombinant IFN-gamma (bovine)	Moderate
Recombinant canine IFN-gamma (canine)	None
Recombinant IFN-gamma (caprine)	Moderate
Recombinant IFN-gamma (chicken)	None
Recombinant IFN-gamma (cynomolgus)	Weak
Recombinant IFN-gamma (dolphin)	Moderate
Recombinant IFN-gamma (equine)	Weak
Recombinant IFN-gamma (feline)	Weak
Recombinant IFN-gamma (human)	Moderate
Recombinant IFN-gamma (mouse)	Weak
Recombinant IFN-gamma (ovine)	Moderate
Recombinant IFN-gamma (rabbit)	Moderate
Recombinant IFN-gamma 1-1 (zebrafish)	None
Recombinant IFN-gamma 1-2 (zebrafish)	None

- > **Limit of detection** calculated as the concentration of IFN-gamma (pig) corresponding to the NSB average plus three standard deviation is  $\leq 0.08$  ng/mL.

## ► Troubleshooting

- **Absorbance values are too low:**
  - one reagent has not been dispensed,
  - incorrect preparation,
  - assay performed before reagents reached room temperature,
  - reading time not long enough.
  
- **High signal and background in all wells:**
  - inefficient washing,
  - overdeveloping (incubation time should be reduced),
  - high ambient temperature.
  
- **High dispersion of duplicates:**
  - poor pipetting technique or irregular plate washing.

These are a few examples of troubleshooting that may occur. If you need further explanation, Bertin Pharma will be happy to assist you. Feel free to contact our technical support staff by phone (+33 (0)139 306 036), fax (+33 (0)139 306 299) or E-mail ([bioreagent@bertinpharma.com](mailto:bioreagent@bertinpharma.com)), and be sure to indicate the batch number of the kit (see outside the box).

Bertin Pharma proposes EIA Training kit #B05005 and EIA workshop upon request. For further information, please contact our Marketing Department by phone (+33 (0)139 306 260) or E-mail ([marketing@bertinpharma.com](mailto:marketing@bertinpharma.com)).



## ► Bibliography

1. Vandenbroeck K., Dijkmans R., van Aerschot A. et al  
Engineering by PCR-based exon amplification of the genomic porcine interferon-gamma DNA for expression in *Escherichia coli*.  
*Biochem Biophys Res Commun.* 1991 Nov 14; 180(3): 1408-15.
2. Billiau A.  
Interferon-gamma: biology and role in pathogenesis.  
*Adv Immunol.* 1996; 62: 61-130.
3. Paludan S.R.  
Interleukin-4 and interferon-gamma: the quintessence of a mutual antagonistic relationship. *Scand J Immunol.* 1998 Nov; 48(5): 459-68.
4. Boehm U., Klamp T., Groot M. et al  
Cellular responses to interferon-gamma.  
*Annu Rev Immunol.* 1997; 15: 749-95
5. Schroder K., Hertzog P.J., Ravasi T. et al  
Interferon-gamma: an overview of signals, mechanisms and functions.  
2004 J. Leukoc. Biol. 75 (2): 163–89.



1	2	3	4	5	6	7	8	9	10	11	12	
A												
B												
C												
D												
E												
F												
G												
H												









Bertin Pharma, over the last decades, has been developing and marketing over 100 biomarker assays, pre-analytical products, kits, antibodies and biochemicals thanks to its innovative work in research and development. Our core areas are orientated to inflammation, oxidative injury, endocrinology, diabetes, obesity, hypertension, neurodegenerative diseases, HIV, prion diseases, pharmacokinetics and metabolism.

Bertin Pharma is active worldwide either with direct sales or through our qualified and trained international distribution network from the United States to Japan.

We are able to provide you with local technical support to use at ease our products.

For further information, please send your request to  
[bioreagent@bertinpharma.com](mailto:bioreagent@bertinpharma.com)



care. innovate. deliver.

Parc d'activités du Pas du Lac - 10 bis avenue Ampère  
F-78180 Montigny-le-Bretonneux - France  
Tel: +33 (0)139 306 036 - Fax: +33 (0)139 306 299  
[bioreagent@bertinpharma.com](mailto:bioreagent@bertinpharma.com) - [www.bertinpharma.com](http://www.bertinpharma.com)

—A **ENIM** Group Company—